


## AQUA:MASTER

The practical solution for irrigating  
the upper and lower gastrointestinal tract

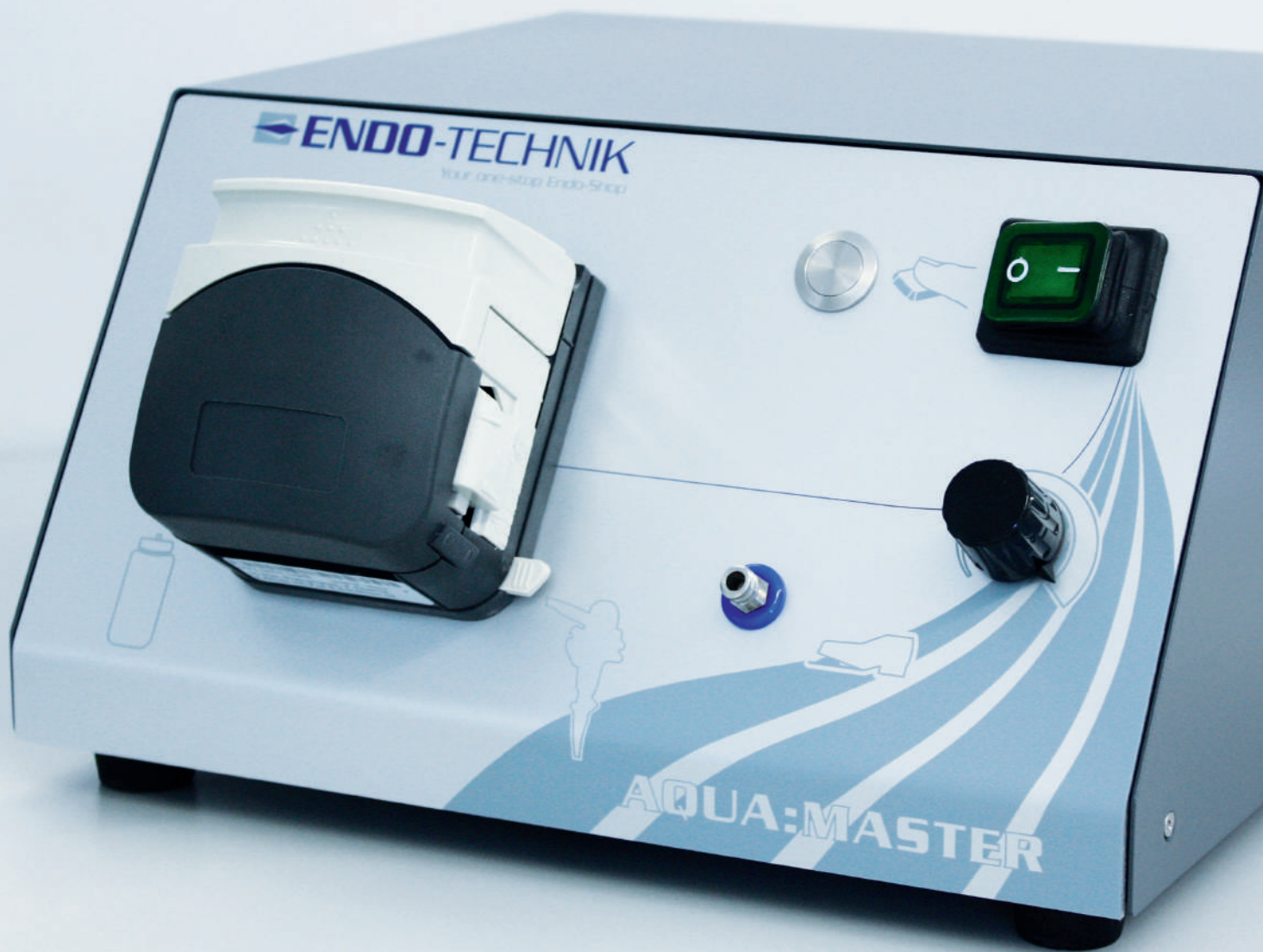

# AQUA:MASTER

## The practical solution for irrigating the upper and lower gastrointestinal tract

The AQUA:MASTER is intended exclusively for irrigating and cleansing the upper and lower gastrointestinal tracts. Direct connection to the water jet or working channel enables the medical practitioner to target intestinal haemorrhaging and remove coagulum and soiling. Our use-once sets of tubes guarantee safe, hygienic application. Disposable sets of tubes with trumpet valves make it possible to perform easy irrigation and extraction even with endoscopes that do not have a water jet channel.



AQUA:MASTER

### Product characteristics

- » Controllable rate of flow: 490 - 1000 ml
- » Simple operation using the hand or foot-operated switch
- » The water tank is suitable for autoclaving
- » Hygienically safe disposable tube system
- » Purging and extraction through Y-valve with no need to disconnect and reconnect

### Technical specifications

Mains voltage:	230 V / 50 Hz
Fuses:	2 x 1.25 A slow-acting
Weight:	5.0 kg
Width:	280 mm
Height:	150 mm
Depth:	300 mm
Certificates:	VDE 0750 (part 1/05.02)

The flow rate is infinitely variable, allowing the volume of fluid used to be adjusted exactly to the application conditions. The low-noise, external peristaltic pump guarantees easy operation and rapid hose changes. Easy operation is ensured by the choice of the hand-operated or foot switch.



AQUA:MASTER rinsing pump



## AQUA:MASTER

14-E400

AQUA:MASTER with foot switch, mains cable, 1 l water bottle, holder for the water bottle, operating instructions, including 10 disposable sets of tubes (14-AD-Y32 suction and pressure tubes) and 14-Y120 Y-adaptor for biopsy channels for all Pentax and Olympus endoscopes with LL biopsy connections, the AQUA:MASTER is supplied with a grey housing from serial number 1437

# AQUA:MASTER

## Accessories for AQUA:MASTER




AQUA:MASTER



Spare part no.: 14-AD-Y32



Spare part no.: 14-AD-Y32-D

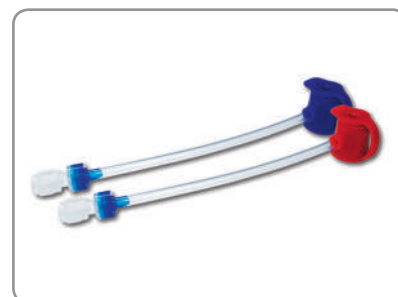

Spare part no.: BIO-VEN-030/ BIO-VEN-050



## AQUA:MASTER ACCESSORIES

14-AD-Y32	Suction and pressure tubes for AQUA:MASTER, 2 x LL-M, overall length 340 cm, of which 70 cm suction tube and 250 cm pressure tube, disposable item
14-AD-Y32-D	Disposable suction and pressure tube for AQUA:MASTER, with pin, 2 x LL-M, overall length 340 cm, of which 70 cm suction tube and 250 cm pressure with flow rate of 470 - 1050 ml, packed in 25s, disposable item
BIO-VEN-030	BPV / biopsy valve with irrigation line for Olympus GI endoscope and Fujifilm series 500, disposable item, blue
BIO-VEN-050	BPV / biopsy valve with irrigation line for all Pentax endoscopes, disposable item, red



Trolleys




Endoscopy Accessories  
General Hospital Equipment



Electronic Accessories  
General Hospital Equipment



Flexible Accessories

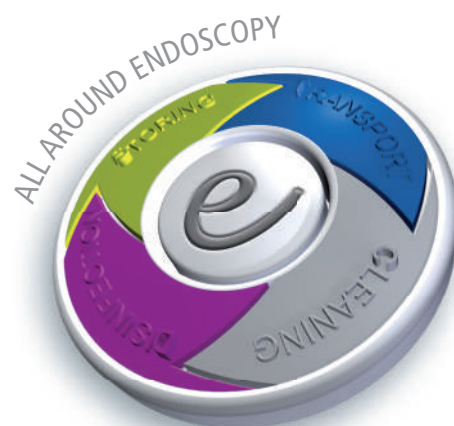

**e-TCDS**

**e**NDOSCOPY

**TRANSPORT**

**CLEANING**

**DISINFECTION**

**STORAGE**

**ACCESSORIES**

**ENDO-TECHNIK W. Griesat GmbH**  
Heiligenstock 32  
42697 Solingen  
Germany

**Tel.:** +49 212 267 43-0  
**Fax:** +49 212 267 43-50

**Email:** [info@endo-technik.de](mailto:info@endo-technik.de)  
**URL:** [www.endo-technik.de](http://www.endo-technik.de)



2017.01

All details, descriptions and images relating to the products are current at the time of going to press these are also subject to change without prior notice. Actual products may differ slightly from images. E&OE.