

Multi-Band-Ligator – 5 Bands

for Esophageal Varices – Single Use – Non Sterile

Multi-Band-Ligator

- User – friendly Mounting
Adaptor – free device fixation onto endoscopes
- Rapid & Ease of Operation of the “Shooting” Mechanism
Facilitates band release in emergency situations
- Click Sound and Visible Indicators
Confirms reliable band release
- Alert Band in Color
Indicates second to last band is placed



Multi-Band-Ligator – 5 Bands

for Esophageal Varices – Single Use – Non Sterile

Medi-Globe's Multi-Band-Ligator features major key advantages for esophageal varices treatment:

➤ Extremely User – friendly Mounting
Adaptor – free device fixation onto endoscopes

➤ Wide Compatibility
Compatible to endoscopes with a tip outer diameter between 8.5 – 11.5 mm

➤ Safe Band Release Procedure
➤ Rapid & Ease of Operation of the “Shooting” Mechanism
Facilitates band releases in emergency situations

➤ Large ergonomic operation wheel
➤ Click Sound and Visible Indicators
Confirms reliable band release

Multi-Band-Ligator – 5 Bands

for Esophageal Varices – Single Use – Non Sterile

➤ Introduction Catheter
Ensures a fast kink-free advancement of the trigger cord through the working channel

➤ Highly Transparent Band Barrel
Perfect varices evaluation

➤ Alert Band in Color
Indicates second to last band on the barrel

➤ Side Port for Irrigation & Flushing

Multi-Band-Ligator – 5 Bands


for Esophageal Varices – Single Use – Non Sterile

Multi-Band-Ligator

Multi-Band Ligator with 5 Bands – Maximum Endoscope Length 150 cm

Article No.	Compatibility	Number of bands	Band inner diameter in released position	Length of the introduction catheter	Min. Working Channel
GBL-01-05-001	Endoscopes with a distal tip diameter between 8.5 – 11.5 mm	5	1.9 mm	150 cm	2.8 mm

Band inner diameter in released status:

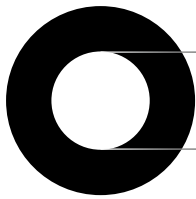

Band inner diameter in released position

Caution: Ligation bands contain natural rubber latex, which may cause allergic reactions.

Distributed by:

Comercializadora Endosistemas Ltda.
Suecia N° 3243, Ñuñoa, Stgo.
Fonos: 2 2209 0262 - 2 2785 5826
ventas@endosistemas.cl

L-GBL-SU-00-ENG Ver. 1.0 (2017-10)