

## Sphincterotomes

- Rotatable Tip  
*Precise Orientation & Positioning*
- 2 & 3 Lumen Versions  
*Contrast Medium Injection while Guide Wire is in Place*
- Monofilament Cutting Wire  
*Accurate & Accurate Cut*
- Version with Semi-insulated Cutting Wire  
*Safe & Controlled Papillary Roof Incision*
- Color Marking with Black Rings  
*Macroscopically Controlled Cannulation Depth*



### Sphincterotomes

Medi-Globe offers a wide range of sphincterotomes in triple, double and single lumen formats. A differentiating feature of Medi-Globe sphincterotomes is the ability to rotate the distal tip from 11 o'clock to 1 o'clock for precise orientation and positioning. Medi-Globe sphincterotomes are available with cutting wires of various lengths and in several tip styles to facilitate biliary and pancreatic duct access in nearly all anatomical requirements.

#### RotaCut®

RotaCut sphincterotomes are equipped with a braided cutting wire.

#### RotaCut® - M

RotaCut - M sphincterotomes are equipped with a monofilament cutting wire.


#### Premier Cut

Premier Cut sphincterotomes are equipped with a "semi-insulated", braided cutting wire.

### S.O.S. Sphincterotomes

A S.O.S. "Speed Option Set" is comprised of a double lumen sphincterotome with braided cutting wire, which is preloaded with a hydrophilic Teflon coated Nitinol guide wire in .035" diameter and a length of 450cm.

#### Sphincterotome Tips

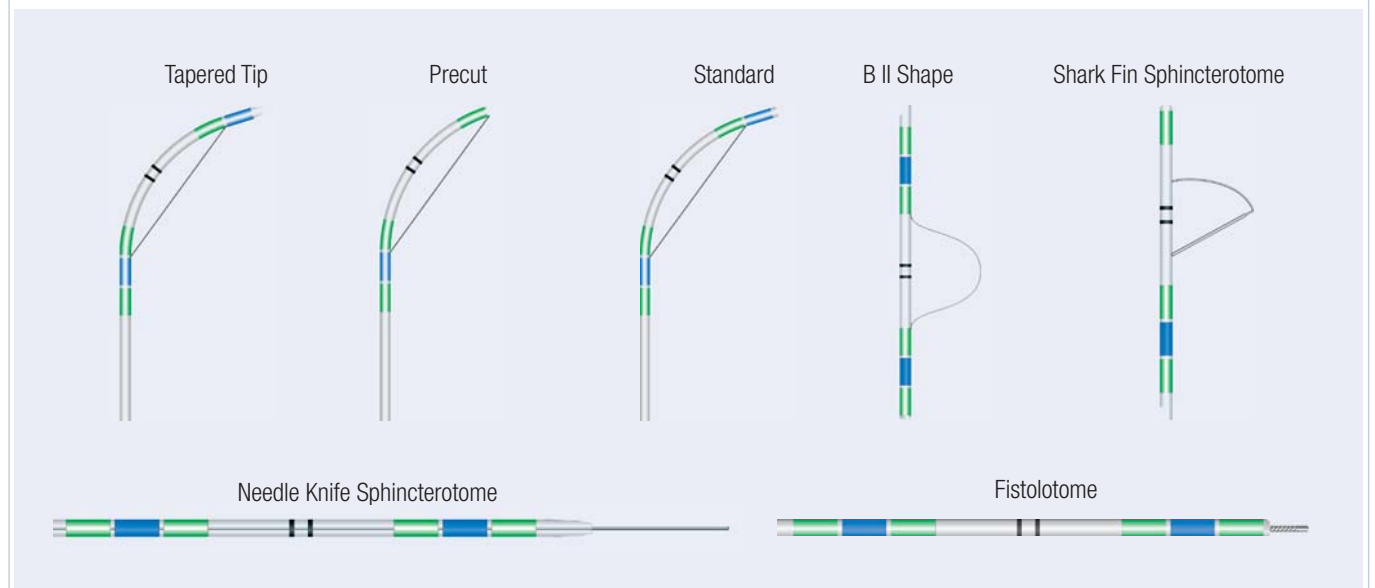

Medi-Globe sphincterotomes are directly compatible with Olympus electro-surgical generators. To verify compatibility with other electro-surgical generators please contact your Medi-Globe sales representative. Please also consult the Medi-Globe "Accessories" leaflet to see a complete list of active cords manufactured to fit most generators.

## 3-Lumen Rotatable Guide Wire Sphincterotomes



Tapered Tip



### RotaCut® – Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-31-20-020	Tapered Tip	5 mm	20 mm	2.55/2.05 mm	.035"
GSP-31-20-025	Tapered Tip	5 mm	25 mm	2.55/2.05 mm	.035"
GSP-31-20-030	Tapered Tip	5 mm	30 mm	2.55/2.05 mm	.035"

### RotaCut® - M – Sphincterotome with Monofilament Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-32-20-020	Tapered Tip	5 mm	20 mm	2.55/2.05 mm	.035"
GSP-32-20-025	Tapered Tip	5 mm	25 mm	2.55/2.05 mm	.035"
GSP-32-20-030	Tapered Tip	5 mm	30 mm	2.55/2.05 mm	.035"

### Premier Cut – Sphincterotome with Braided and Semi-insulated Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-27-20-030	Tapered Tip	5 mm	20 mm	2.55/2.05 mm	.035"

## 1-Lumen Rotatable Sphincterotomes



Needle Knife



### Needle Knife Sphincterotomes

Catalog No.	Needle Projection	Outer Sheath Diameter
GSP-11-17-005	0 - 5 mm	1.7 mm
GSP-11-17-010	0 -10 mm	1.7 mm
GSP-11-17-015	0 -15 mm	1.7 mm

## 2-Lumen Rotatable Guide Wire Sphincterotomes



Precut



### RotaCut® – Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-22-20-020	Precut Tip	-	20 mm	2.05/2.55 mm	.035"
GSP-22-20-030	Precut Tip	-	30 mm	2.05/2.55 mm	.035"
GSP-22-17-020	Precut Tip	-	20 mm	1.7 mm	.021"
GSP-22-17-030	Precut Tip	-	30 mm	1.7 mm	.021"

B II



### RotaCut® – Billroth II Sphincterotome with Braided Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Max. Accepted Guide Wire
GSP-34-20-020	Billroth II	-	Variable	2.05/2.55 mm	.035"

Tapered Tip



### Premier Cut – Sphincterotome with Braided and Semi-insulated Cutting Wire

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Accepted Guide Wire
GSP-25-20-020	Tapered Tip	5 mm	20 mm	2.05/2.55 mm	.035"
GSP-25-20-030	Tapered Tip	5 mm	30 mm	2.05/2.55 mm	.035"
GSP-26-20-020	Tapered Tip	15 mm (long nose)	20 mm	2.05/2.55 mm	.035"
GSP-26-20-030	Tapered Tip	15 mm (long nose)	30 mm	2.05/2.55 mm	.035"

Tapered Tip



### S.O.S. Sphincterotome – "Speed Option Set" with Pre-loaded Hydrophilically Teflon Coated Nitinol Guide Wire in .035" Diameter and 450 cm Length

Catalog No.	Description	Tip Length	Cutting Wire Length	Outer Sheath Diameter	Accepted Guide Wire
GSP-21-20-031	RotaCut®, Tapered Tip	5 mm	30 mm	2.55 / 2.05 mm	Pre-loaded
GSP-26-20-031	Premier Cut, Semi-insulated Cutting Wire, Tapered Tip	5 mm	30 mm	2.55 / 2.05 mm	Pre-loaded

#### Distributed by:

Comercializadora Endosistemas Ltda.  
Suecia N° 3243, Ñuñoa, Stgo.  
Fono: 22090262 - F/Fax: 27855826  
ventas@endosistemas.cl

#### Manufactured By:

Medi-Globe GmbH Medi-Globe-Straße 1-5 83101 Achenmühle, Germany  
☎ +49 80 32 / 973- 379 📠 +49 80 32 / 973-339 ✉ sales@medi-globe.de  
[www.medi-globe.de](http://www.medi-globe.de)