

AUTOMATIC ENDOSCOPE REPROCESSOR

ENDO CLEAN 1000




Endoscope is delivered under most reliable environment with strong cleaning capability and sterilization ability

- Regular model of ENDO CLEAN 1000, with its multiple 3 channel nozzle injection type , provides strong cleaning capability for endoscope and also prevents mixing with sterilization liquid by removing as much moisture as possible remaining in endoscope by deposition type of draining.
- It's also possible to set time freely taking average 7 to 15 minutes
- Sterilization is available for all types of endoscope and 1 set endoscope is possible for sterilization.
- Material for Basin is chemical-resistance indicating excellent durability while it has another benefit making it possible to install endoscope without a mounting device.




Providing effective sterilization capability with washing liquid and alcohol added

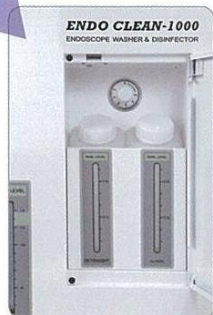

- Maintaining concentration of sterilization liquid for certain period by minimizing dilution rate
- A function to prevent draining of sterilization liquid
- Easy to determine sterilization liquid replacement cycle with counter display of sterilization liquid
- Using parts with strong chemical resistance against Peracetic Acid and PAA
- Automatic alcohol injection at inside channel

Strong multiple 3 channel 360° rotating nozzle type

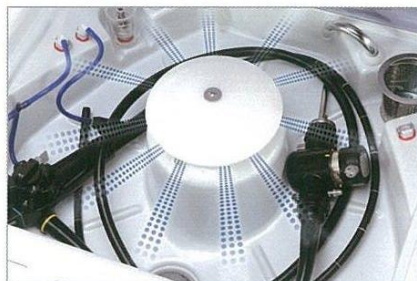


Program mode selection type which is easy to be operated
[washable and sterilized with minimum time]



Equipped with strong air gun

[Drying external skin and internal hollow rapidly ↑] [Preventing nozzle clogging ↑]

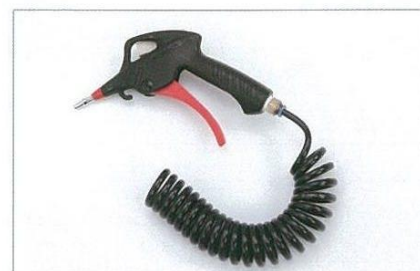

Providing OLYMPUS, PENTAX, FUJINON leak jack adaptor



Automatic Process



Applicable Scopes	Flexible endoscopes
Visual Leakage Detection	Bubble detection during immersion
Number of Reprocessor Endoscope	for 1 scope
Washing(Detergent) Method	
Exterior surfaces	Multi channel jetting system rotated in a 360° arc
Channel interiors	Fluid flushing cleaning
Accessories	Fluid flushing cleaning
Disinfection Method	
Exterior surfaces	Disinfectant solution immersion
Channel interiors	Disinfectant solution flushing and filling
Accessories	Disinfectant solution immersion
Washing Time	10sec ~ 3min
Detergent Time	5sec
Disinfection Time	1min ~ 1hr
Alcohol flushing Time	5sec
Disinfectant Tank (Built-in heater) Temperature setting	15°C ~ 50°C
Tank capacity	Disinfectant 15L Detergent 1L Alcohol 1L
Demension (mm)	550(W) X 585(D) X 950(H)
Weight	75Kg
Power	220-240VAC 50/60Hz

